

# Invacare® | **Birdie® and Birdie® Compact**



5 year warranty\*!



**think**  
safe patient handling



NEW

NEW feature -  
180kg (28 stone)  
weight limit

## Transportability and patient comfort go hand-in-hand

With maximum space for the user, both the *Invacare Birdie* and the *Birdie Compact* offer a truly comfortable lift and transfer to or from beds, chairs or even the floor. The lifters are designed to ensure that folding and unfolding can be carried out easily and without the need for tools. In addition, the lifter can be dismantled into two parts if needed – also without tools.



**Yes, you can.®**



**Total user comfort**

The high boom maximizes space for the client even when the hoist is in the highest position. This gives an ergonomic lift and transfer as the client can be easily rotated in the hoist. This makes placing the client into a bed/wheelchair much easier for the carer. The optimization of space creates a reduction in risk of injury as the client's knees will not be close to the actuator of the hoist.

**Ease of handling: simple to fold, unfold and dismantle without tools**

Handling the lifter has been made as smooth and easy as possible as the hoist can easily be folded and unfolded without the need for tools. The **Birdie** takes up minimal space when stored and is easy to push and transport. If required, the mast and boom can be easily detached from the base (no tools required).

**TÜV-certification for maximum security**

For total security, the hoists come with a snap-lock system, thus minimizing the risk of accident. The lifters are CE-certified in line with the Medical Device Directive 93/42/ECC. To ensure

the hoists comply with EN ISO 10535 (the European standard for hoists) the hoist has been rigorously tested both by Invacare and the German TÜV institute, an accredited external body.

**Features and Options**



**Long boom – high user comfort**

The user can easily be rotated and positioned. No risk of injuries to the knees.



**Easy to fold and dismantle**

Ease of transport and storage. No tools required when folding / unfolding or when dismantling the hoist.



**Snap lock**

Optimal security to avoid risk of accidents.



**Emergency stop and emergency lowering**

Manual or electric emergency lowering (compact version features only electric lowering). Available with fixed or detachable battery. Off board battery charger is available.

**Features and Options**



**Click on spreader bar**  
Easy to change the spreader bar.



**Mounted spreader bar**  
Correctly and safe positioning of the spreader bar.



**Manual leg spread**  
Simply step on the pedal to adjust the leg spread.



**Handle for leg spread**  
Handle can easily be attached to operate the manual leg spread.



**Electrical leg spread**  
Available for *Invacare Birdie* (option).



**2 point quick release spreader bars**  
Easy to attach and detach. Available widths: 350 mm, 450 mm, 550 mm.




**Push bar and hand control**  
Superior ergonomics for carers.



**Castors**  
100 mm or 75 mm castors available. (braking castors are shown)



**Scale**  
If required, a scale is available approved. According to the medical standard EN 45501, certified and sealed in order to be used to medicate.

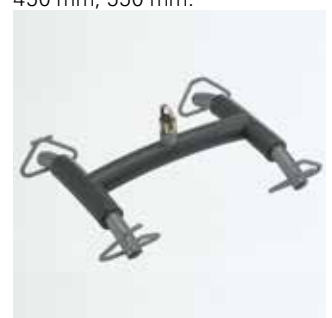

**4 point quick release spreader bars**  
Easy to attach and detach. Available widths: 450 mm and 550 mm.

## Technical data

|                       |        |       |       |                 |       |
|-----------------------|--------|-------|-------|-----------------|-------|
|                       |        |       |       |                 |       |
| <b>Birdie</b>         | 124 cm | 64 cm | 42 kg | 44,5 - 171 cm   | 10 cm |
| <b>Birdie Compact</b> | 109 cm | 52 cm | 36 kg | 52,5 - 161,5 cm | 10 cm |

|                       |         |        |           |        |        |
|-----------------------|---------|--------|-----------|--------|--------|
|                       |         |        |           |        |        |
| <b>Birdie</b>         | 101 cm  | 180 kg | 40 cycles | 2,9 Ah | 140 cm |
| <b>Birdie Compact</b> | 84,5 cm | 150 kg | 40 cycles | 2,9 Ah | 107 cm |

All measures are taken with 75 mm castors and 2 point spreaderbar.  
For 100 mm castors add 15 mm for height measures and 20 mm for width.

## Frame colour



Light grey /  
Dark grey

## \*Warranty:

Actuator: 5 years  
Frame: 3 years  
Battery: 1 year

Invacare Limited  
Pencoed Technology Park - Pencoed -  
Bridgend CF35 5AQ - United Kingdom  
Tel: +44 1 656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail:  
[ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Ireland Ltd  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS - County  
Dublin - Ireland  
Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [Ireland@invacare.com](mailto:Ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)

© 2014 Invacare International Sàrl  
All rights reserved.

All information quoted is believed to be  
correct at time of print.  
Invacare reserves the right to alter  
product specifications without prior  
consultation.

Birdie - UK - 11/2014 - 1504431